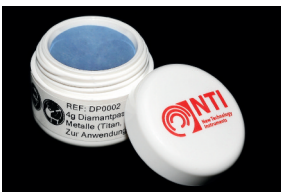












Application Chart • Brushes for the dental laboratory

for an efficient application of polishing brushes



| | Natural bristles | Natural bristles | Composite polishing brush | Natural bristles | Natural bristles | Felt polishers FiPol | Cotton buff | Abrasive buffs SoftPol | |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---|
| opt. Speed | 15.000 rpm | 15.000 rpm | 6.000 rpm | 10.000 rpm | 10.000 rpm | 10.000 rpm | 3.000 rpm | 10.000 rpm | |
|  <p>REF DP0002 Contents 4g</p> <ul style="list-style-type: none"> ▲ Pre Polishing ■ Refined Finish ● High-Shine Polishing <p>* Use with Diamond polishing paste</p> |  <p>P1259 P1261 P1268 soft</p> |  <p>P1260 P1262 P1267</p> |  <p>P1270</p> |  <p>P1271 hard</p> |  <p>P1272 extra-hard</p> |  <p>P4123 P4124 Density extra-hard</p> |  <p>P1269</p> |  <p>P4201 coarse</p>  <p>P4202 medium</p>  <p>P4203 fine</p> | |
| Denture acrylics | | | | | | | | | |
| Partial and complete dentures | ▲ | ▲ | | ▲ | ▲ | | ● | ▲ | |
| Composites | | | | | | | | | |
| Veneers | ▲ | ▲ | ▲ | ▲ | | | ● | | |
| Dental ceramics* | | | | | | | | | |
| Facets, veneers, inlays, bridges | ▲ | ▲ | | ▲ | ▲ | ■ | ● | | |
| Precious alloys | | | | | | | | | |
| Inlays, crown, bridges | ▲ | ▲ | | ▲ | ▲ | | ● | | |
| Semi-precious alloys | | | | | | | | | |
| Inlays, crowns, bridges | ▲ | ▲ | | ▲ | ▲ | | ● | | |
| Non-precious alloys | | | | | | | | | |
| Crowns, bridges, Milling technique* | ▲ | ▲ | | ▲ | ▲ | ■ | ● | | |
| Metal frame work* | | | | | | | | | |
| | ▲ | ▲ | | ▲ | ▲ | | ● | | |
| Titanium | | | | | | | | | |
| Crowns, bridges, metal frame work | ▲ | ▲ | | ▲ | ▲ | | ● | | |
| Soft relining material | | | | | | | | | ▲ |
| Vacuum-formed foils | | | | | | | | | ▲ |